

Is There a Simple and Easy Way to Detect Singlet Oxygen? Comparison of Methods for Detecting Singlet Oxygen and Application to Measure Scavenging Activity of Various Compounds

Tokuko Takajo, Kazunori Anzai*

Faculty of Pharmaceutical Sciences, Nihon Pharmaceutical University; 10281 Komuro, Ina-machi, Kita-Adachi-gun, Saitama 362-0806, Japan

*Correspondence should be addressed to Kazunori Anzai; anzai@nichiyaku.ac.jp

Received date: October 24, 2020, **Accepted date:** December 22, 2020

Copyright: © 2020 Takajo T, et al. This is an open-access article distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited.

Abstract

Reactive oxygen species (ROS) are known to exert both beneficial and harmful effects in the human body. Singlet oxygen (1O_2), is highly reactive and considered as one of the ROS, although it is not a radical molecule. 1O_2 reacts with many kinds of biological components such as lipids, proteins, and nucleic acids. 1O_2 is short lived because it reacts with biomolecules and collisions with water molecules rapidly causing the return of 1O_2 to the ground state, and is therefore not easy to quantify. Antioxidants, such as β -carotene, lycopene and astaxanthin, are quenchers of 1O_2 . To compare the 1O_2 quenching activity of new test compounds with these typical quenchers, it is desirable to measure 1O_2 specifically, speedily, and easily. In this paper, various methods for the measurement of 1O_2 , such as UV-Vis spectroscopy, fluorescence spectroscopy, and electron spin resonance, were compared for the purpose of evaluating the 1O_2 quenching activity of various test compounds.

Keywords: Singlet oxygen (1O_2), 1O_2 detection probe, Quenching activity, UV-Vis spectroscopy, Fluorescence spectroscopy, Near-infrared spectroscopy, Electron spin resonance

Introduction

Reactive oxygen species (ROS) are known to exert both beneficial and harmful effects in the human body. The beneficial aspects of ROS are, for example, roles in signaling cascades and immunological defense. In contrast, overproduction of ROS has been reported to be involved in various diseases such as cancer, cardiovascular disease, and diabetes. To prevent the harmful effects of ROS, organisms have antioxidant enzymes and various antioxidative compounds. Various antioxidative compounds are present in foods, beverages, agricultural products, plants, as well as in the human body [1]. Measurement of antioxidant activity is important for developing new drugs, health foods, and nutritional supplements.

Atmospheric oxygen (triplet oxygen, 3O_2) is a stable form of oxygen molecule. The other form, singlet oxygen

(1O_2), is highly reactive and considered as one of the ROS, although it is not a radical molecule. 1O_2 reacts with many kinds of biological components such as lipids, proteins, and nucleic acids [2-4]. 1O_2 is short lived because it reacts with biomolecules and collisions with water molecules rapidly causing the return of 1O_2 to the ground state. Thus, it is not easy to quantify 1O_2 . 1O_2 is thought to be a cause of some disorders such as photo-aging, skin damage, and erythropoietic porphyria [5]. On the other hand, cell death and tissue destruction induced by 1O_2 are used for photodynamic therapy (PDT), a promising cancer treatment [6-8].

It is known that antioxidants, such as β -carotene, lycopene, and astaxanthin, are quenchers of 1O_2 . In order to compare the 1O_2 quenching activity of new test compounds with these typical quenchers, it is desirable to specifically, speedily, and easily measure 1O_2 . In addition, since light and heat may decompose the structure and

affect the quenching activity of the test compound, it is not ideal to use light or heat for $^1\text{O}_2$ measurement. In this paper, we summarize and compare various methods for the measurement of $^1\text{O}_2$ for the purpose of evaluating the $^1\text{O}_2$ quenching activity of various test compounds. Artificial $^1\text{O}_2$ generation is usually performed either by photosensitization of photosensitizers, such as rose bengal and porphyrins, or thermal decomposition of 4-methyl-1,4-etheno-2,3-benzodioxin-1(4*H*)-propanoic acid (EP).

UV-Vis Spectroscopy

UV-Vis absorbance measurement is simple, speedy, and the most commonly used method. Anthracene and benzofuran are known as compounds that selectively scavenge $^1\text{O}_2$. The absorbance of these aromatic molecules decreases when they are oxidized by $^1\text{O}_2$, and the amount of $^1\text{O}_2$ can be estimated semi-quantitatively from the decrease in absorbance. UV-Vis spectroscopy is somewhat less sensitive than fluorescence spectroscopy. However, UV-Vis spectroscopy is less prone to artifacts, since the absorption spectrum of a compound with high absorption coefficient cannot be dominated by small amounts of contaminants. One of the most widely used $^1\text{O}_2$ probes for UV-Vis spectroscopy is 1,3-diphenylisobenzofuran (DPBF, Figure 1), which is commercially available at a low cost. DPBF reacts with $^1\text{O}_2$ in an essentially diffusion-limited reaction rate and thus has very high sensitivity. It traps up to 50% of $^1\text{O}_2$ produced in alcohol/water or

micellar solution. In addition, replacing H_2O with D_2O (a $^1\text{O}_2$ lifetime extender), is an option when only small amounts of $^1\text{O}_2$ are produced [9]. However, DPBF is not specific to $^1\text{O}_2$ and readily reacts with other ROS such as hydroxyl radical, alkoxy radical ($\text{RO}\cdot$), and H_2O_2 [10-12]. Furthermore, if the test compound is colored, the color affects the absorbance measurement with DPBF.

Fluorescence Spectroscopy

Comparing the quantitative measurements of absorbance and fluorescence, the detection limit of fluorescence is generally three orders of magnitude lower. To detect $^1\text{O}_2$, a fluorescent sensor called Singlet Oxygen Sensor (SOSG, Figure 1) [13] has been used widely, ranging from material studies to medical applications, although the probe is expensive. SOSG exhibits weak blue fluorescence with excitation peaks at 372 nm and 393 nm and emission peaks at 395 nm and 416 nm. Upon reaction with $^1\text{O}_2$, endoperoxide, the immediate product of SOSG, exhibits green fluorescence with excitation and emission peaks at 504 nm and 525 nm, respectively. Therefore, SOSG can readily detect $^1\text{O}_2$. However, it was demonstrated that SOSG can generate $^1\text{O}_2$ under exposure to UV (355 nm) or visible light (532 nm) in the absence of photosensitizers [14]. It was confirmed that the irradiation of visible light induces SOSG transitions to a triplet excited state, which converts $^3\text{O}_2$ to $^1\text{O}_2$ [15]. Therefore, a positive error can occur with respect to the amount of $^1\text{O}_2$ measured by SOSG. In

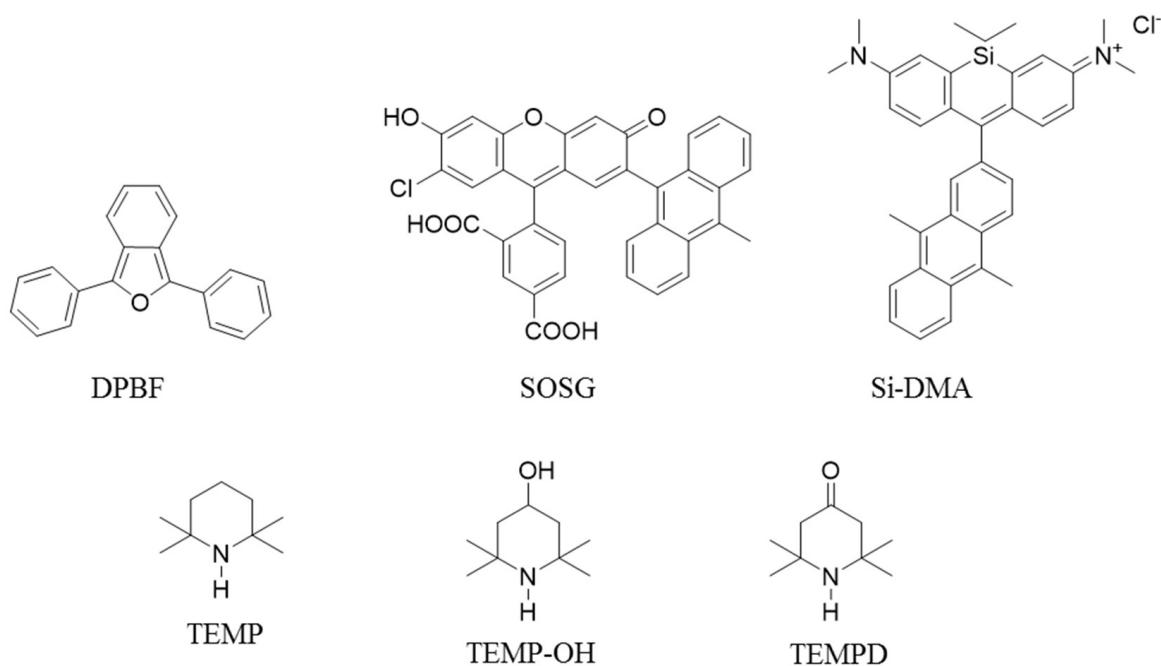


Figure 1: Chemical structure of detector probes.

addition, since many photosensitizers used to generate $^1\text{O}_2$ have a pronounced absorption band in the visible region, if the absorption overlaps with the emission peak of the fluorophores of probes, the signal produced by the reaction with $^1\text{O}_2$ will be diminished. Photosensitizers, such as rose bengal and eosin Y, exhibit an intense absorption band at 525 nm, which coincides with the emission peak of SOSG endoperoxide [16]. Incorrect use of SOSG may lead to misinterpretation of the data on the yield of $^1\text{O}_2$ production by photosensitization reactions. For biological application, a new type of fluorescence probe for $^1\text{O}_2$ (Si-DMA, Figure 1) was developed [17]. Si-DMA has an emission peak at 680 nm and is permeable through cell membranes. It is now commercially available and can be used for imaging of $^1\text{O}_2$ produced in mitochondria and by PDT.

Near-infrared Spectroscopy

$^1\text{O}_2$ can be directly measured by detecting weak light emission in the near infrared region. Phosphorescence at 1270 nm emitted by $^1\text{O}_2$ can be detected specifically at a low quantum yield [18]. By using an improved near-infrared-sensitive photomultiplier tube, time-resolved $^1\text{O}_2$ luminescence was measured in various solutions of aluminum tetrasulphonated phthalocyanine and Photofrin [19]. Time-resolved analysis showed a strongly reduced lifetime of $^1\text{O}_2$ and an increased lifetime of triplet state of the photosensitizer in the intracellular component [19]. This method is often used for imaging in cells and animals. However, the method requires specialized instruments because of the low quantum yield.

Singlet Oxygen Absorption Capacity

Recently, a new assay method to quantify $^1\text{O}_2$ absorption capacity (SOAC) was developed [20-23]. The SOAC values were measured in ethanol/chloroform/ D_2O (50:50:1, v/v/v) solution at 35°C using a UV-Vis spectrophotometer equipped with a six-channel cell positioner and a temperature control unit. The reaction of antioxidants with $^1\text{O}_2$ was measured using the competition reaction method, where EP was used for $^1\text{O}_2$ generation and DPBF was used as an absorbance probe. In the SOAC method, $^1\text{O}_2$ quenching activity is evaluated kinetically, and the activity of various phenolic antioxidants can be compared. The method is useful for evaluating the $^1\text{O}_2$ quenching activity of both lipophilic and hydrophilic antioxidants. The method allows for measurement of activity that differs in five orders of magnitude (different rate constants of $10^{10} \text{ M}^{-1}\text{s}^{-1}$ for lycopene to $10^5 \text{ M}^{-1}\text{s}^{-1}$ for ferulic acid). However, in the SOAC assay, decrease in the absorbance is observed in 120 to 180 min, during which the test compound may become photodegraded.

Electron Spin Resonance

Electron spin resonance (ESR) is a non-destructive

method that detects unpaired electrons and is specific to free radicals. Since $^1\text{O}_2$ is a ROS but is not a free radical, it cannot be detected directly by ESR. A $^1\text{O}_2$ detector probe, 2,2,6,6-tetramethylpiperidine (TEMP, Figure 1), was introduced in 1976 [24]. TEMP reacts with $^1\text{O}_2$ to produce a stable nitroxide radical, 2,2,6,6-tetramethylpiperidine-*N*-oxyl (TEMPO), which can be quantified by ESR spectroscopy. Other sterically hindered secondary amines, such as 2,2,6,6-tetramethyl-4-piperidinol (TEMP-OH, Figure 1) and 2,2,6,6-tetramethyl-4-piperidone (TEMPD, Figure 1), are also used as $^1\text{O}_2$ detector probes [25-27]. Photoirradiation of rose bengal is often used as a $^1\text{O}_2$ generation method. When $^1\text{O}_2$ generated in the rose bengal system is detected by ESR using a hindered secondary amine probe, the signal of the nitroxide radical derived from the secondary amine probe is measured. However, rose bengal reportedly generates ROS other than $^1\text{O}_2$, such as super oxide anion radical and hydroxyl radical [26,28], which may affect the ESR signal. Indeed, when visible light was irradiated to the mixed solution of TEMP-OH and rose bengal, obvious splitting was observed in the ESR spectrum. This observation suggested the overlap of 4-oxo-2,2,6,6-tetramethylpiperidine-*N*-oxyl (TEMPO) and 4-hydroxy-2,2,6,6-tetramethylpiperidine-*N*-oxyl (TEMPOL) [25]. It was reported that the exposure of TEMPOL to hydroxyl radical led to the appearance of a new triplet ESR signal attributable to TEMPO, along with a reduction in the intensity of TEMPOL signals [29]. The following is a possible explanation for the ESR signal splitting observed in the mixed solution of rose bengal and TEMP-OH. TEMP-OH reacts with $^1\text{O}_2$ to form TEMPOL, which is oxidized to TEMPO by hydroxyl radicals generated from rose bengal. On the other hand, when EP was used for $^1\text{O}_2$ generation and TEMP-OH was used as a detector probe, no splitting of the three-line ESR spectrum was obtained, although the signal intensity was very small [25]. When the mixture solution of EP and TEMP-OH was reacted at 35°C for 40 min, the ESR signal intensity was smaller than that obtained with irradiation of rose bengal with visible light for 3 min, supporting that $^1\text{O}_2$ production from EP is smaller than that from rose bengal. For the detection by ESR, if EP is used as a $^1\text{O}_2$ generator, it is recommended to add D_2O in the system to increase the signal intensity. Although the amount of $^1\text{O}_2$ generated from EP is small, it is possible to test the $^1\text{O}_2$ scavenging activity of test compounds with a reaction time of 30~40 min in the presence of D_2O using TEMP-OH as a detector probe.

Conclusion

Our opinions with respect to the pros and cons of the detection methods for $^1\text{O}_2$ are summarized in Table 1. When measuring the scavenging activity of various compounds to $^1\text{O}_2$, it is desirable to avoid light irradiation, because the scavenging activity may decrease due to the decomposition of the test compounds. Since the ESR method does not

Methods	Pros	Cons	Probes
UV-Vis spectroscopy	Easy	Low sensitivity Disturbance by color	DPBF
Fluorescent spectroscopy	High sensitivity	Disturbance by color	SOSG Si-DMA
Near-infrared spectroscopy	Direct measurement High specificity	Low sensitivity Needs special equipment	
Singlet oxygen absorption capacity	Multi sample measurement	Photodegradation of sample Disturbance by color	DPBF
Electron spin resonance	High specificity No disturbance by color	Needs ESR spectrometer	TEMP TEMP-OH TEMPD

Table 1: Comparison of the methods for detecting $^1\text{O}_2$.

require light irradiation, it is suitable for measuring the $^1\text{O}_2$ scavenging activity of various compounds. In addition, EP is better for the $^1\text{O}_2$ generation system than the photosensitization method, because EP does not require light irradiation. In the present paper, we only focused on the $^1\text{O}_2$ measurement system. However, development of a new $^1\text{O}_2$ generation system may also contribute to improving the methods for measuring the $^1\text{O}_2$ scavenging activity of various compounds. Because the methodologies described so far may detect other ROS together, it is necessary to develop specific probes capable of detecting $^1\text{O}_2$ only, or specific reagents or systems that generate only $^1\text{O}_2$ with moderate amounts. Finding compounds or plants with new $^1\text{O}_2$ scavenging capacity by measuring scavenging activity correctly and easily should be useful to prevent or treat photoaging, skin damage, and erythropoietic porphyria caused by $^1\text{O}_2$.

References

- Mangels AR, Holden JM, Beecher GR, Forman MR, Lanza E. Carotenoid content of fruits and vegetables: an evaluation of analytic data. *Journal of the American Dietetic Association.* 1993 Mar 1;93(3):284-96.
- Di Mascio P, Martinez GR, Miyamoto S, Ronsein GE, Medeiros MH, Cadet J. Singlet molecular oxygen reactions with nucleic acids, lipids, and proteins. *Chemical Reviews.* 2019 Feb 5;119(3):2043-86.
- Davies MJ, Truscott RJ. Photo-oxidation of proteins and its role in cataractogenesis. *Journal of Photochemistry and Photobiology B: Biology.* 2001 Oct 1;63(1-3):114-25.
- Davies MJ. Singlet oxygen-mediated damage to proteins and its consequences. *Biochemical and Biophysical Research Communications.* 2003 Jun 6;305(3):761-70.
- Kon T, Tanigawa T, Hayamizu K, Shen M, Tsuji T, Naito Y, et al. Singlet oxygen quenching activity of human serum. *Redox Report.* 2004 Dec 1;9(6):325-30.
- Weishaupt KR, Gomer CJ, Dougherty TJ. Identification of singlet oxygen as the cytotoxic agent in photo-inactivation of a murine tumor. *Cancer Research.* 1976 Jul 1;36(7 Part 1):2326-9.
- Dolmans DE, Fukumura D, Jain RK. Photodynamic therapy for cancer. *Nature Reviews Cancer.* 2003 May;3(5):380-7.
- Henderson BW, Waldow SM, Mang TS, Potter WR, Malone PB, Dougherty TJ. Tumor destruction and kinetics of tumor cell death in two experimental mouse tumors following photodynamic therapy. *Cancer Research.* 1985 Feb 1;45(2):572-6.
- Entradas T, Waldron S, Volk M. The detection sensitivity of commonly used singlet oxygen probes in aqueous environments. *Journal of Photochemistry and Photobiology B: Biology.* 2020 Mar 1;204:111787.
- Carloni P, Damiani E, Greci L, Stipa P, Tanfani F, Tartaglini E, Wozniak M. On the use of 1, 3-diphenylisobenzofuran (DPBF). Reactions with carbon and oxygen centered radicals in model and natural

systems. *Research on Chemical Intermediates.* 1993 Jan 1;19(5):395-405.

11. Ohyashiki T, Nunomura M, Katoh T. Detection of superoxide anion radical in phospholipid liposomal membrane by fluorescence quenching method using 1, 3-diphenylisobenzofuran. *Biochimica et Biophysica Acta (BBA)-Biomembranes.* 1999 Sep 21;1421(1):131-9.

12. Żamojć K, Zdrowowicz M, Rudnicki-Velasquez PB, Krzyminiński K, Zaborowski B, Niedziałkowski P, et al. The development of 1, 3-diphenylisobenzofuran as a highly selective probe for the detection and quantitative determination of hydrogen peroxide. *Free Radical Research.* 2017 Jan 2;51(1):38-46.

13. Flors C, Fryer MJ, Waring J, Reeder B, Bechtold U, Mullineaux PM, et al. Imaging the production of singlet oxygen in vivo using a new fluorescent sensor, Singlet Oxygen Sensor Green®. *Journal of Experimental Botany.* 2006 May 1;57(8):1725-34.

14. Ragàs X, Jiménez-Banzo A, Sánchez-García D, Batllori X, Nonell S. Singlet oxygen photosensitisation by the fluorescent probe Singlet Oxygen Sensor Green®. *Chemical Communications.* 2009(20):2920-2.

15. Kim S, Fujitsuka M, Majima T. Photochemistry of singlet oxygen sensor green. *The Journal of Physical Chemistry B.* 2013 Nov 14;117(45):13985-92.

16. DeRosa MC, Crutchley RJ. Photosensitized singlet oxygen and its applications. *Coordination Chemistry Reviews.* 2002 Nov 1;233:351-71.

17. Kim S, Tachikawa T, Fujitsuka M, Majima T. Far-red fluorescence probe for monitoring singlet oxygen during photodynamic therapy. *Journal of the American Chemical Society.* 2014 Aug 20;136(33):11707-15.

18. Ogilby PR. Singlet oxygen: there is indeed something new under the sun. *Chemical Society Reviews.* 2010;39(8):3181-209.

19. Niedre M, Patterson MS, Wilson BC. Direct Near-infrared Luminescence Detection of Singlet Oxygen Generated by Photodynamic Therapy in Cells In Vitro and Tissues In Vivo¶. *Photochemistry and Photobiology.* 2002 Apr;75(4):382-91.

20. Ouchi A, Aizawa K, Iwasaki Y, Inakuma T, Terao J, Nagaoka SI, et al. Kinetic study of the quenching reaction of singlet oxygen by carotenoids and food extracts in solution. Development of a singlet oxygen absorption capacity (SOAC) assay method. *Journal of Agricultural and Food Chemistry.* 2010 Sep 22;58(18):9967-78.

21. Aizawa K, Iwasaki Y, Ouchi A, Inakuma T, Nagaoka SI, Terao J, et al. Development of singlet oxygen absorption capacity (SOAC) assay method. 2. Measurements of the SOAC values for carotenoids and food extracts. *Journal of Agricultural and Food Chemistry.* 2011 Apr 27;59(8):3717-29.

22. Mukai K, Ouchi A, Takahashi S, Aizawa K, Inakuma T, Terao J, et al. Development of singlet oxygen absorption capacity (SOAC) assay method. 3. Measurements of the SOAC values for phenolic antioxidants. *Journal of Agricultural and Food Chemistry.* 2012 Aug 15;60(32):7905-16.

23. Iwasaki Y, Takahashi S, Aizawa K, Mukai K. Development of singlet oxygen absorption capacity (SOAC) assay method. 4. Measurements of the SOAC values for vegetable and fruit extracts. *Bioscience, Biotechnology, and Biochemistry.* 2015 Feb 1;79(2):280-91.

24. Lion Y, Delmelle M, Van de Vorst A. New method of detecting singlet oxygen production. *Nature.* 1976 Sep;263(5576):442-3.

25. Takajo T, Kurihara Y, Iwase K, Miyake D, Tsuchida K, Anzai K. Basic Investigations of Singlet Oxygen Detection Systems with ESR for the Measurement of Singlet Oxygen Quenching Activities. *Chemical and Pharmaceutical Bulletin.* 2020 Feb 1;68(2):150-4.

26. Nakamura K, Ishiyama K, Ikai H, Kanno T, Sasaki K, Niwano Y, Kohno M. Reevaluation of analytical methods for photogenerated singlet oxygen. *Journal of Clinical Biochemistry and Nutrition.* 2011;49:87-95.

27. Nardi G, Manet I, Monti S, Miranda MA, Lhiaubet-Vallet V. Scope and limitations of the TEMPO/EPR method for singlet oxygen detection: the misleading role of electron transfer. *Free Radical Biology and Medicine.* 2014 Dec 1;77:64-70.

28. Criado S, Allevi C, Ceballos C, García NA. Visible-light promoted degradation of the commercial antioxidants butylated hydroxyanisole (BHA) and butylated hydroxytoluene (BHT): a kinetic study. *Redox Report.* 2007 Dec 1;12(6):282-8.

29. Saito K, Takeshita K, Ueda JI, Ozawa T. Two reaction sites of a spin label, TEMPOL (4-hydroxy-2, 2, 6, 6-tetramethylpiperidine-N-oxyl), with hydroxyl radical. *Journal of Pharmaceutical Sciences.* 2003 Feb 1;92(2):275-80.